


Terumo *Glidecath*® 4 Fr

Unprecedented smoothness and access to small peripheral vessels

- Large lumen (0.038" wire compatible) and small profile (4 Fr) is ideal for:
 - use as a guiding catheter for microcatheters, such as Progreat™ Microcatheters
 - diagnostic procedures that require high flow rates
- Flexible construction for excellent trackability and navigation of even the most tortuous anatomies

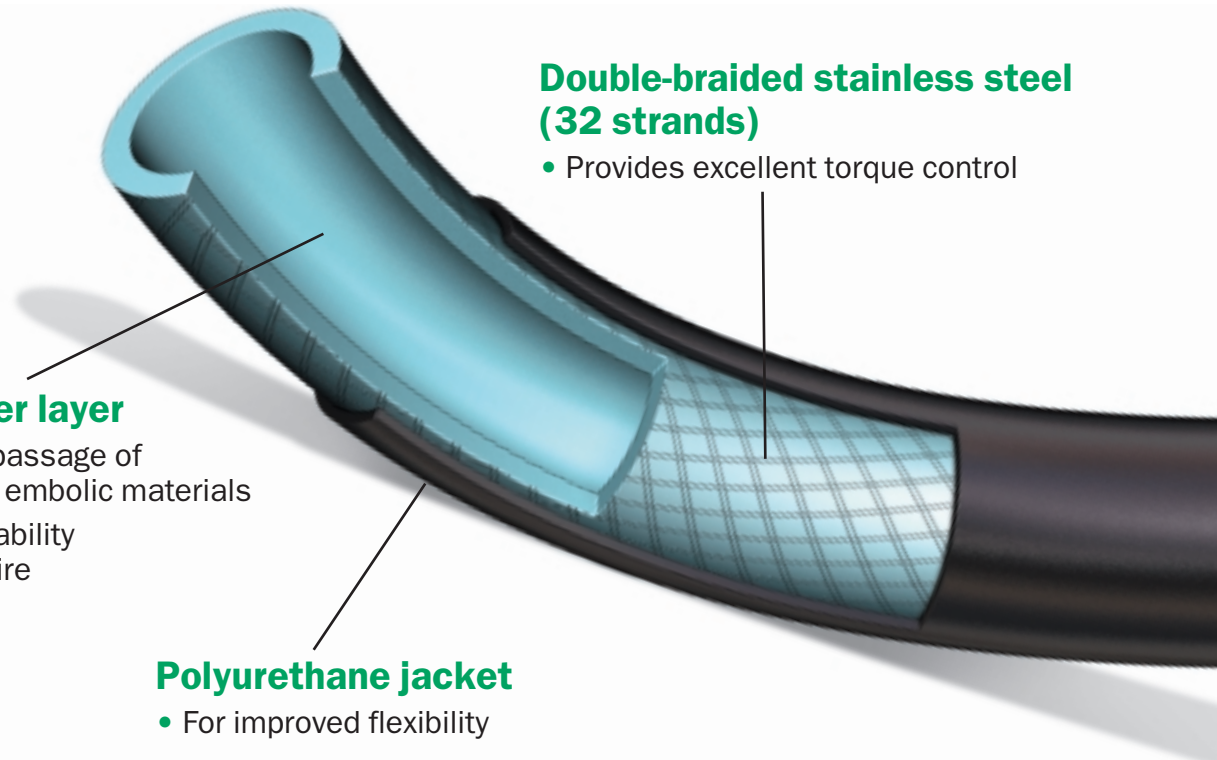

Double-braided stainless steel (32 strands)


- Provides excellent torque control

Nylon-rich inner layer

- Allows smooth passage of therapeutic and embolic materials
- Enhances trackability over the guidewire

Polyurethane jacket

- For improved flexibility

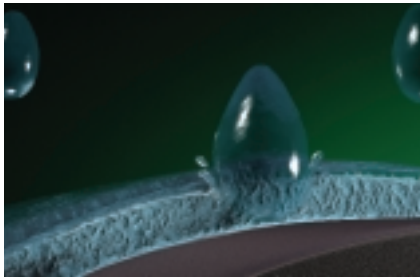

Terumo Glide Technology™ hydrophilic coating on distal 15 cm

- Unique hydrophilic coating has been imitated but never equaled
- Enables smooth navigation and more efficient handling












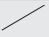






Glidecath and Glidewire® together — unsurpassed performance

Hydrophilic Coated Catheters

Terumo *Glidecath*® 4 Fr

Hydrophilic Coated Catheters

0.038" guidewire compatible and packed 5 per box.

Product Code	Catheter Size (Fr)	Total Length (cm)	Flow Rate* (mL/sec)	Hydrophilic Coating Length (cm)	Tip Shape	
CG412	4	65	15	15	Straight	
CG415	4	65	15	15	Angle	
CG408	4	65	15	15	C1	
CG409	4	65	15	15	C2	
CG411	4	65	15	15	J-Tip	
CG413	4	100	13	15	Straight	
CG416	4	100	13	15	Angle	
CG410	4	100	13	15	C2	
CG418	4	100	13	15	MP	
CG401	4	100	13	15	Sim 1	
CG402	4	100	13	15	Sim 2	
CG403	4	100	13	15	Sim 3	
CG404	4	100	13	15	H1	
CG405	4	100	13	15	JB1	
CG406	4	100	13	15	JB2	
CG407	4	100	13	15	Mani	
CG414	4	120	11	15	Straight	
CG417	4	120	11	15	Angle	

*All measurements taken at 750 psi and 37°C using Renografin® 76.

For more information, call 800.862.4143 to speak to a Customer Care Representative, or visit www.terumomedical.com/interventional
 Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2006 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. Lit. # TIS011 04/06