


Terumo *Glidecath*® XP 5 Fr

Durability that goes the distance

- Maintains its integrity during long procedures
- Large lumen (0.038" wire compatible) is ideal for:
 - use as a guiding catheter for microcatheters, such as Progreat™ Microcatheters
 - diagnostic procedures that require high flow rates

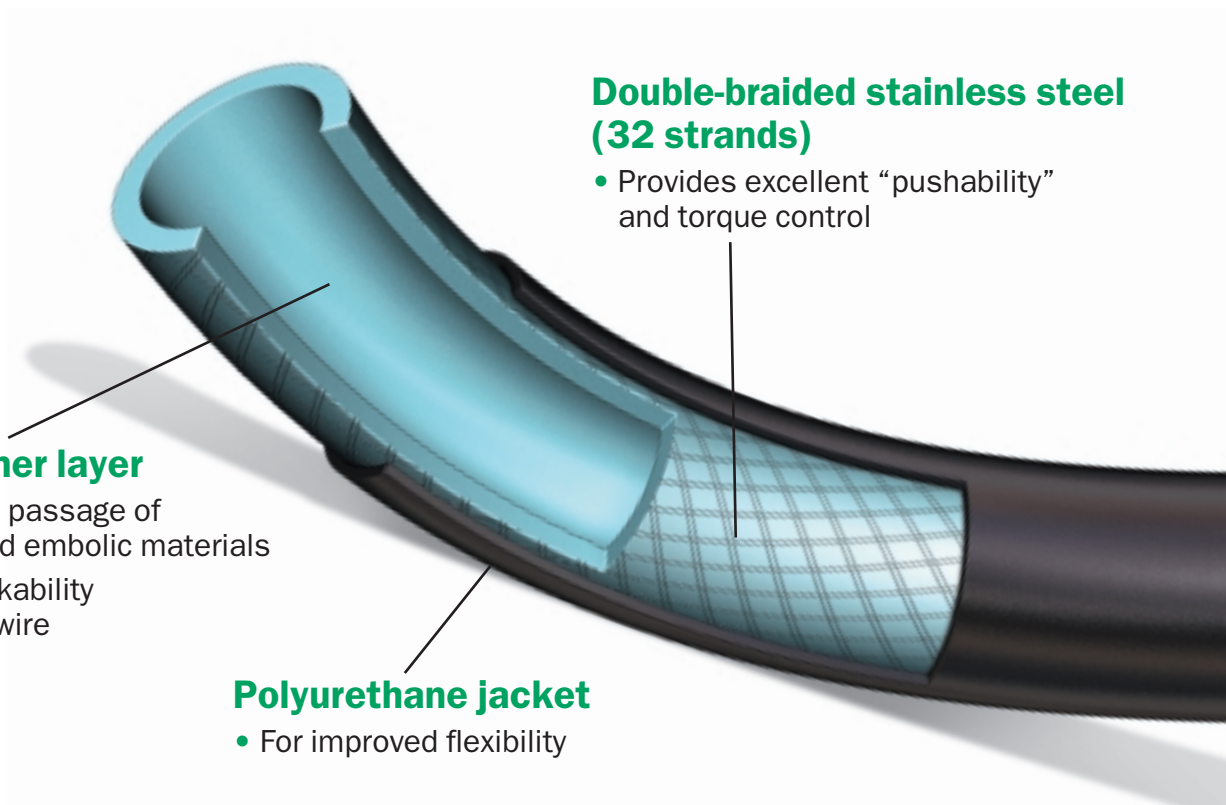

Double-braided stainless steel (32 strands)


- Provides excellent “pushability” and torque control

Nylon-rich inner layer

- Allows smooth passage of therapeutic and embolic materials
- Enhances trackability over the guidewire

Polyurethane jacket

- For improved flexibility

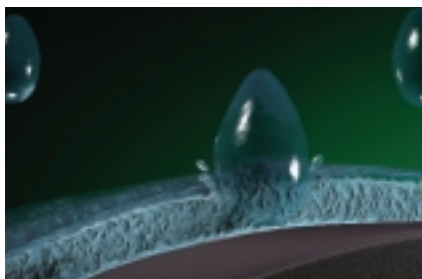

Terumo Glide Technology™ hydrophilic coating on distal 15 cm

- Unique hydrophilic coating has been imitated but never equaled
- Enables smooth navigation and more efficient handling









Glidecath and Glidewire® together — unsurpassed performance

Hydrophilic Coated Catheters

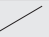Terumo *Glidecath*® XP 5 Fr

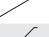Hydrophilic Coated Catheters

0.038" guidewire compatible and packed 5 per box.

Product Code	Catheter Size (Fr)	Flow Rate* (mL/sec)	Total Length (cm)	Hydrophilic Coating Length (cm)	Tip Shape	
RF*ZM7510GM	5	24	65	15	Straight	
RF*ZV9510GM	5	24	65	15	Angle	
RF*ZB5510GM	5	24	65	15	C2	
RF*ZM75110M	5	21	100	15	Straight	
RF*ZV95110M	5	21	100	15	Angle	
RF*WA15110M	5	21	100	15	S1	
RF*WA25110M	5	21	100	15	S2	
RF*WB15110M	5	21	100	15	H1	

*All measurements taken at 1000 psi and 37°C using Renografin® 76.

For more information, call 800.862.4143 to speak to a Customer Care Representative, or visit www.terumomedical.com/interventional
 Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2006 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. Lit. # TIS013 04/06