


# Terumo *Glidecath*® XP 5 Fr

## Durability that goes the distance

- Maintains its integrity during long procedures
- Large lumen (0.038" wire compatible) is ideal for:
  - use as a guiding catheter for microcatheters, such as Progreat™ Microcatheters
  - diagnostic procedures that require high flow rates

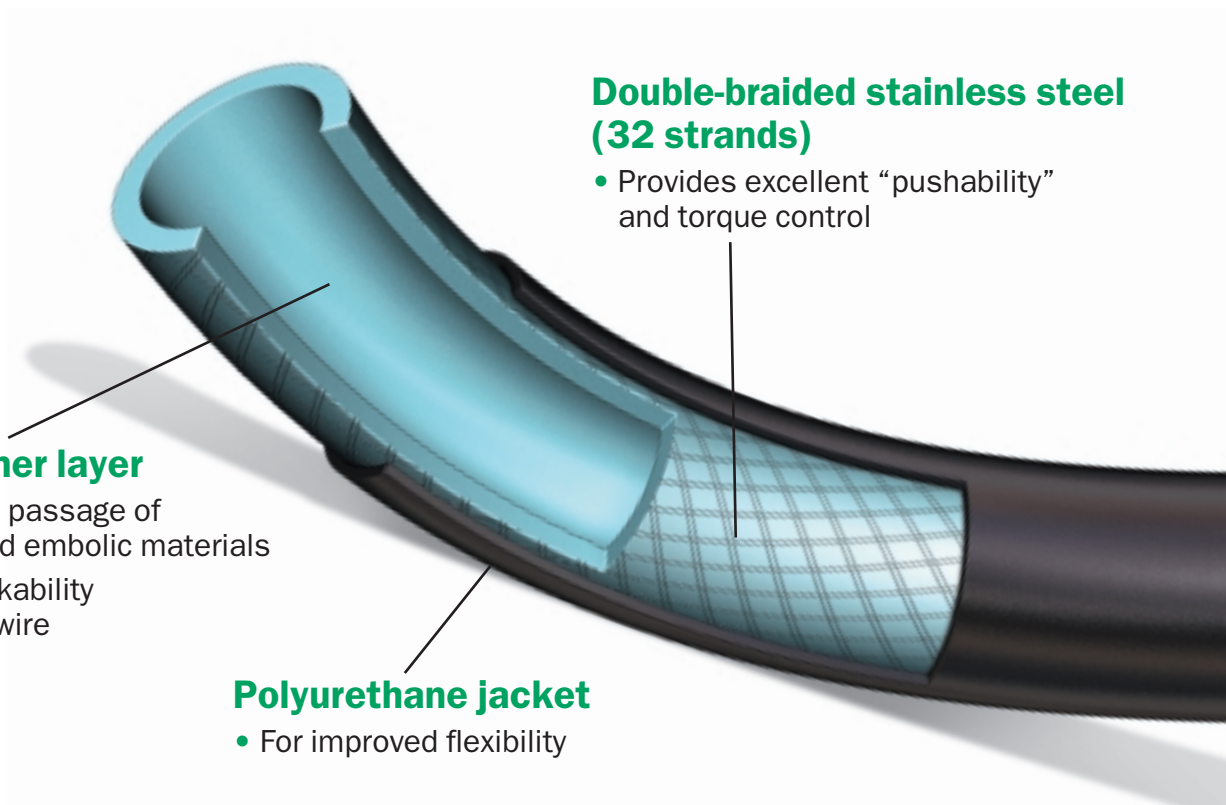

### Nylon-rich inner layer


- Allows smooth passage of therapeutic and embolic materials
- Enhances trackability over the guidewire

### Double-braided stainless steel (32 strands)

- Provides excellent “pushability” and torque control

### Polyurethane jacket

- For improved flexibility

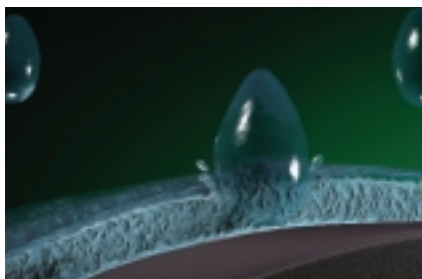

### Terumo Glide Technology™ hydrophilic coating on distal 15 cm

- Unique hydrophilic coating has been imitated but never equaled
- Enables smooth navigation and more efficient handling









**Glidecath and Glidewire® together — unsurpassed performance**

## Hydrophilic Coated Catheters
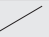
# Terumo *Glidecath*® XP 5 Fr
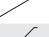
### Hydrophilic Coated Catheters

0.038" guidewire compatible and packed 5 per box.

Product Code	Catheter Size (Fr)	Flow Rate* (mL/sec)	Total Length (cm)	Hydrophilic Coating Length (cm)	Tip Shape	
RF*ZM7510GM	5	24	65	15	Straight	
RF*ZV9510GM	5	24	65	15	Angle	
RF*ZB5510GM	5	24	65	15	C2	
RF*ZM75110M	5	21	100	15	Straight	
RF*ZV95110M	5	21	100	15	Angle	
RF*WA15110M	5	21	100	15	S1	
RF*WA25110M	5	21	100	15	S2	
RF*WB15110M	5	21	100	15	H1	

\*All measurements taken at 1000 psi and 37°C using Renografin® 76.

For more information, call 800.862.4143 to speak to a Customer Care Representative, or visit [www.terumomedical.com/interventional](http://www.terumomedical.com/interventional)  
 Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2006 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. Lit. # TIS013 04/06