

Terumo *Glidewire*® GT

Go where you need to go. See what you need to see.

- 25 cm distal taper minimizes vessel trauma and improves selectivity of vessels
- Wide selection of tip angles, including the unique double angle, to navigate even the most acutely angled vessels

Terumo Glide Technology™ hydrophilic coating

- Unmatched lubricity consistently gets you where you need to go

Terumo nitinol alloy core

- Resists kinking for easier, faster navigation through even the most tortuous vessels
- Restores shape to manufactured state for multiple use in a single procedure

2 cm gold coil at distal tip


- Provides enhanced radiopacity of tip for optimal visibility

Polyurethane jacket with tungsten


- Provides enhanced radiopacity of entire Terumo Glidewire

Core-to-tip design

- Provides optimal torque control for superior navigation

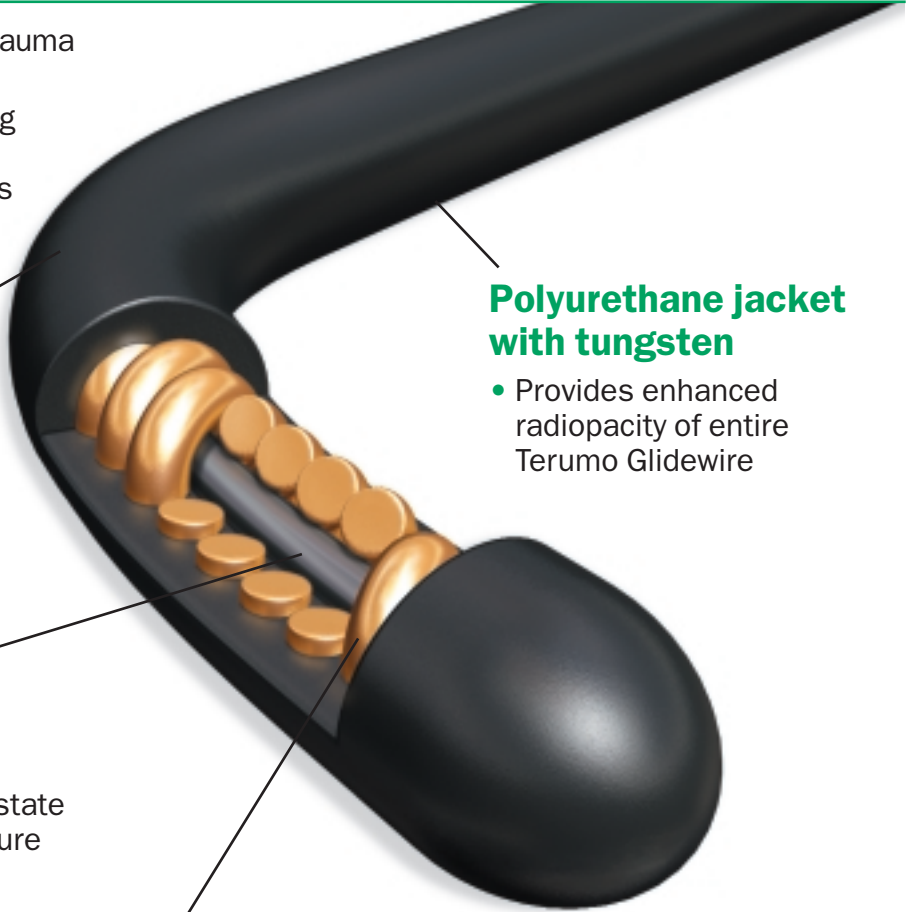

Consistent, reliable Terumo Glidewire performance that gets you there

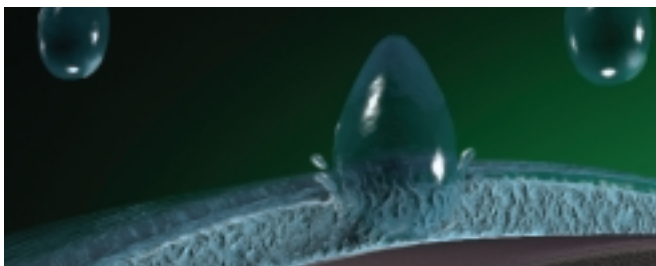

Terumo Glide Technology hydrophilic coating





- Unique hydrophilic coating has been imitated but never equaled
- Consistent performance gets you there faster

Small Vessel Guidewires (0.018")

Terumo *Glidewire*® GT

Super-Selective Hydrophilic Coated Guidewires

Packed 1 per box and include a torque device.

Product Code	Diameter (in)	Total Length (cm)	Flexible Tip Type	Grind Length (cm)	Radiopaque Length (mm)	Tip Shape
RG*GA1818SA	0.018	180	Standard	25	20	45° Angle 
RG*GE1818SA	0.018	180	Standard	25	20	90° Angle 
RG*GW1818SA	0.018	180	Standard	25	20	90°/60° Double Angle [†] 
RG*ES1818SA	0.018	180	Standard	25	20	Straight/Shapeable 

[†]Distal/proximal.

For more information, call 800.862.4143 to speak to a Customer Care Representative, or visit www.terumomedical.com/interventional
Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2006 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. Lit. #TIS005 04/06