

Terumo *Glidewire*[®] Gold/*Glidewire*[®] Gold Neuro

The Terumo *Glidewire* you know and trust with
“Gold Standard” visibility for the most demanding procedures

Glidewire Gold

- Choose from 3 cm or 5 cm distal taper

Glidewire Gold Neuro

- 8 cm distal taper

Terumo Glide Technology[™] hydrophilic coating

- Unmatched lubricity consistently gets you where you need to go

Core-to-tip design

- Provides optimal torque control for superior navigation

Polyurethane jacket with tungsten

- Provides enhanced radiopacity of entire Terumo *Glidewire*


Terumo nitinol alloy core

- Resists kinking for easier, faster navigation through even the most tortuous vessels
- Restores shape to manufactured state for multiple use in a single procedure

2 mm gold band at distal tip

- Provides the enhanced radiopacity you need in smaller vessels

*Consistent, reliable Terumo *Glidewire* performance that gets you there*

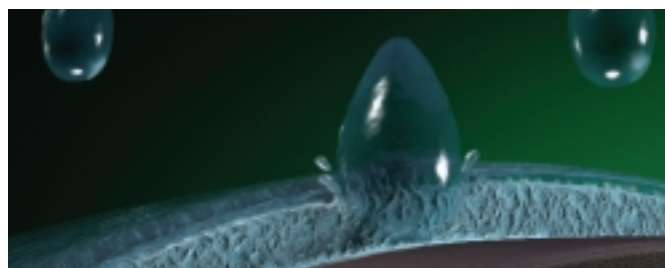

Terumo Glide Technology hydrophilic coating

- Unique hydrophilic coating has been imitated but never equaled
- Consistent performance gets you there faster

Small Vessel Guidewires (0.018")

Terumo *Glidewire*® Gold

Hydrophilic Coated Guidewires

Packed 1 per box and include a torque device.

Product Code	Diameter (in)	Total Length (cm)	Flexible Tip Length (cm)	Radiopaque Length (mm)	Tip Shape
GM1811	0.018	180	3	2	45° Angle
GM1813	0.018	180	3	2	70° Angle
GM1812	0.018	180	5	2	45° Angle
GM1814	0.018	180	5	2	70° Angle

Neuro Guidewires (0.011"-0.018")

Terumo *Glidewire*® Gold Neuro

Hydrophilic Coated Guidewires for Neuro Applications

Packed 1 per box and include a torque device.

Product Code	Diameter (in)	Total Length (cm)	Flexible Tip Length (cm)	Radiopaque Length (mm)	Tip Shape
GM1001	0.011	180	8	2	45° Angle
GM1002	0.011	180	8	2	70° Angle
GM1401	0.014	180	8	2	45° Angle
GM1402	0.014	180	8	2	70° Angle
GM1601	0.016	180	8	2	45° Angle
GM1602	0.016	180	8	2	70° Angle
GM1801	0.018	180	8	2	45° Angle
GM1802	0.018	180	8	2	70° Angle

For more information, call 800.862.4143 to speak to a Customer Care Representative, or visit www.terumomedical.com/interventional
Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2006 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. Lit. #TIS004 04/06