

Pinnacle® R/O II

Radiopaque Marker Introducer Sheaths (4 Fr - 9 Fr)

Packed 10 per box and includes dilator.

	Product Code	Length (cm)	Inner Diameter (Fr)	Dilator Protruding Length (cm)	Mini Guidewire Diameter (in)
Pinnacle R/O II without Mini Guidewire	RSB401	6	4	2.5	-
	RSB501	6	5	2.5	-
	RSB601	6	6	2.5	-
	RSB701	6	7	2.5	-
	RSB801	6	8	2.5	-
	RSB901	6	9	2.5	-
	RSB402	10	4	2.5	-
	RSB502	10	5	2.5	-
	RSB602	10	6	2.5	-
	RSB702	10	7	2.5	-
	RSB802	10	8	2.5	-
	RSB902	10	9	2.5	-
	RSB403	25	4	2.5	-
	RSB503	25	5	2.5	-
	RSB603	25	6	2.5	-
Pinnacle R/O II with Mini Guidewire	RSB703	25	7	2.5	-
	RSB803	25	8	2.5	-
	RSB903	25	9	2.5	-
	RSB411	6	4	2.5	0.035
	RSB511	6	5	2.5	0.038
	RSB611	6	6	2.5	0.038
	RSB711	6	7	2.5	0.038
	RSB811	6	8	2.5	0.038
	RSB911	6	9	2.5	0.038
	RSB412	10	4	2.5	0.035
	RSB512	10	5	2.5	0.038
	RSB612	10	6	2.5	0.038
RSB712	10	7	2.5	0.038	
RSB812	10	8	2.5	0.038	
RSB912	10	9	2.5	0.038	

Maximum guidewire size: 0.035" for 4 Fr; 0.038" for all others.

Terumo — A Global Leader in Medical Devices

Introducers

Pinnacle® Introducer Sheaths
Pinnacle® Destination® Guiding Sheaths
Pinnacle® R/O II Radiopaque Marker Introducer Sheaths
Pinnacle® R/O II HiFlo Introducer Sheaths
Pinnacle TIF Tip™ Introducer Sheaths
Glidesheath™ Access Kits
Glidesheath™ Introducer Sheaths
GlideAccess™ Micro Accessing Systems

Guidewires

Glidewire® Hydrophilic Coated Guidewires
Glidewire® GT Hydrophilic Coated Guidewires
Glidewire® Gold Hydrophilic Coated Guidewires
Runthrough® NS Coronary Guidewires

Angiographic Catheters

Glidecath® Hydrophilic Coated Catheters
Glidecath® XP Hydrophilic Coated Catheters

Coronary Diagnostic Catheters

Optitorque™ Diagnostic Catheters

Micro-Guide Catheters

FineCross™ MG Coronary Micro-Guide Catheters

Microcatheters

Progreat™ Hydrophilic Coated Microcatheters
Progreat™ Coaxial Microcatheter Systems

Radial Compression Devices

TR Band™ Radial Compression Devices

Interventional Products

AZUR™ Framing Coil Embolization Systems
AZUR™ Peripheral HydroCoil® Embolization Systems
Bead Block™ Compressible Microspheres

Accessories

Torque Device
Obturator

For more information, call 800.862.4143 to speak to an Inside Sales Customer Care Representative, or visit www.terumo.com
Terumo Interventional Systems · 2101 Cottontail Lane · Somerset, NJ 08873 · Fax: 800.411.5870

©2009 Terumo Medical Corporation. All rights reserved. HydroCoil is a registered trademark of MicroVenton Inc. Bead Block is a trademark of Biocompatibles UK Ltd. All other brand names are trademarks or registered trademarks of Terumo. Lit.TIS021 02/09

This is TERUMO Territory.™


Pinnacle® R/O II

Radiopaque Marker Introducer Sheaths

Available in 4 Fr - 9 Fr sizes to meet all of your procedural needs

It's easy to see how it's easy to see

- Platinum/iridium band shows precise location of distal tip

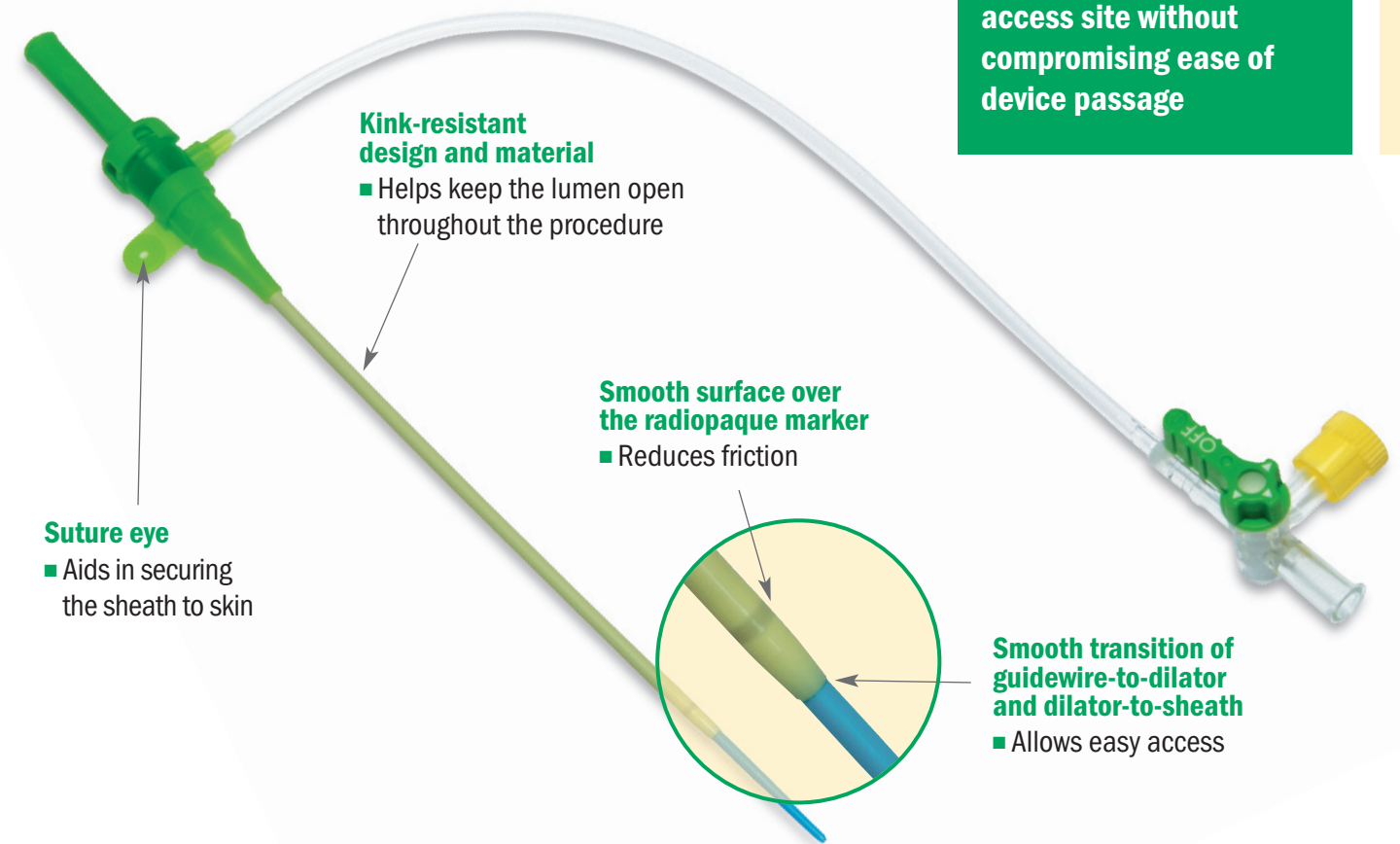

Small outside diameter and large inside diameter minimizes the size of the access site without compromising ease of device passage

When access makes all the difference, This is Terumo Territory.™

The recognized Pinnacle sheath has performance advantages you have come to trust



- Unique Terumo Cross-Cut Valve maintains uncompromised hemostasis
- "Snap-on" dilator lock prevents backout during insertion
- Hub has suture eye for securing sheath to skin and is color-coded for accurate sizing