


Glidecath[®]

Hydrophilic Coated Catheter

ENGAGE AND DELIVER

Expand your access and delivery options
for therapeutic delivery

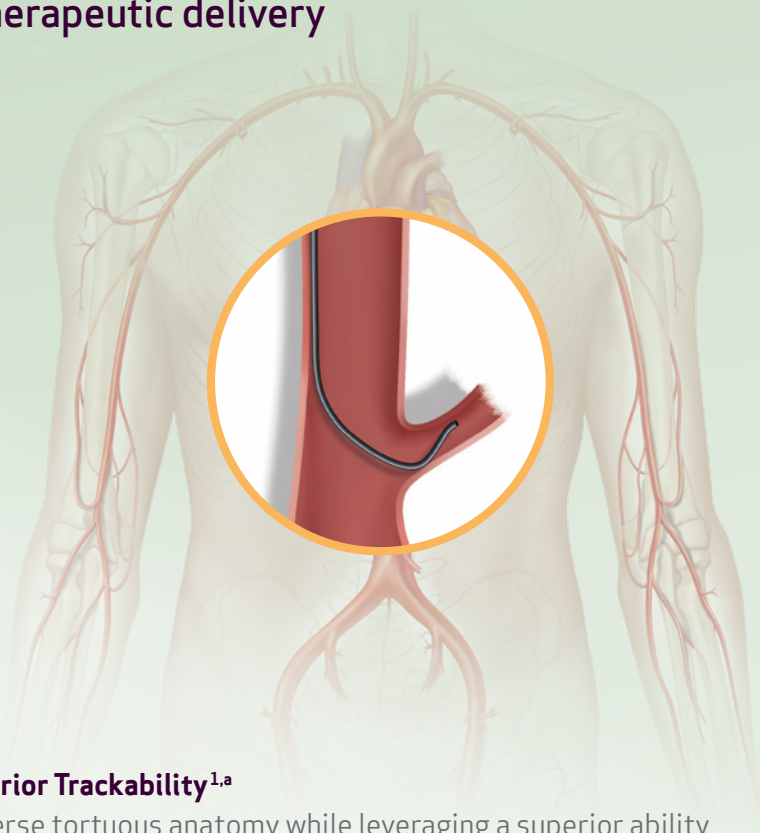

- **Superior Trackability^{1,a}**
Traverse tortuous anatomy while leveraging a superior ability to deliver devices through a nylon rich inner layer
- **Superior Torqueability^{1,a}**
Engage acute vessels with superior torque control with a double-braided design
- **Excellent Lubricity^{1,a}**
Navigate difficult vasculature with the excellent lubricity provided by Terumo Glide Technology

^a When compared to leading competitors

PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you achieve **better outcomes for more patients.**



GLIDECATH® Hydrophilic Coated Catheter MG1 and MG2

PRODUCT CODE	FRENCH SIZE	LENGTH (CM)	TIP SHAPE	HYDROPHILIC COATING (CM)
CG440	4Fr	125 cm	MG1	Distal 15 cm
CG441	4Fr	135 cm	MG2	Distal 40 cm
CG442	4Fr	150 cm	MG2	Distal 40 cm
RF*ZRV511CM	5Fr XP	125 cm	MG1	Distal 15 cm
RF*ZUF511DJM	5Fr XP	135 cm	MG2	Distal 40 cm
RF*ZUF5115JM	5Fr XP	150 cm	MG2	Distal 40 cm



GLIDECATH® tip shapes specifically designed for radial access provides you with the performance required for visceral and pelvic interventions.

GLIDECATH® Hydrophilic Coated Catheter MG1 and MG2 Specifications

	4Fr	5Fr XP
Effective Lengths	125 cm, 135 cm, 150 cm	125 cm, 135 cm, 150 cm
Construction	Double-Braided Stainless Steel	Double-Braided Stainless Steel
Inner Diameter - Distal Tip	0.0405" / 1.03 mm	0.0472" / 1.20mm
Inner Diameter - Proximal Shaft	0.0413" / 1.05 mm	0.0480" / 1.22 mm
Guidewire Compatibility	0.038"	0.038"
Maximum Injection Pressure	750 psi	1000 psi
Hydrophilic Coating	Distal 15 cm and 40 cm	Distal 15 cm and 40 cm

FIND OUT MORE



US: 800.862.4143
CA: 833.883.7866



US: terumo.com
CA: terumocanada.ca

Reference:
1. Data on file.

RX ONLY. The advertisement is directed to physicians only, and not to consumers. Refer to product labels and packaging insert for complete warnings, precautions, potential complications, and instructions for use. Products may not have regulatory approval in all countries. Please contact your local sales representative if you have questions about the availability of products in your area.

©2021 Terumo Medical Corporation and Terumo Medical Canada Inc. All rights reserved. The third-party trademarks used herein are trademarks of their respective owners. All other brand names are trademarks or registered trademarks of Terumo. PM-04314

