

# ***FineCross® M3***

Coronary Micro-Guide Catheter

# ***Runthrough® NS Izanai™***

Coronary Guidewires

**UNIQUELY  
ENGINEERED**  
FOR SIMPLE TO  
COMPLEX PCI



 **TERUMO**  
INTERVENTIONAL  
SYSTEMS



# PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you apply effective solutions and achieve **better outcomes for more patients.**

 **TERUMO**  
INTERVENTIONAL  
SYSTEMS




# FineCross® M3

Coronary Micro-Guide Catheter

## OPTIMAL GUIDEWIRE SUPPORT AND SUPERIOR CROSSABILITY<sup>1</sup>

FINECROSS® M3 Coronary Micro-Guide Catheter is designed to provide the ultimate balance of properties for access and crossing of complex lesions.


- 1.7 Fr tapered distal tip is designed for improved crossability and guidewire handling<sup>1</sup>
- Next Generation Hydrophilic Coating, with its unmatched and most durable lubricity, provides excellent trackability and crossability<sup>1</sup>
- 15 cm floppy distal segment allows for improved trackability
- Tapered stainless steel braid construction is designed for optimal guidewire support

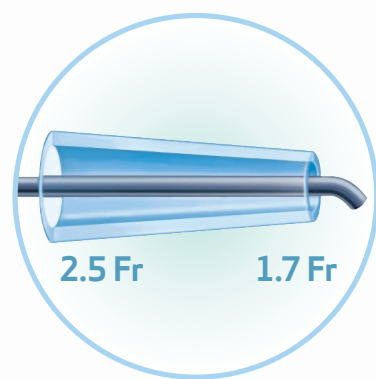

### Tapered Distal Tip

Tapered ID and OD along the length of the catheter for unsurpassed crossability<sup>2</sup>

Outer proximal diameter: 2.5 Fr  
Outer distal diameter: 1.7 Fr

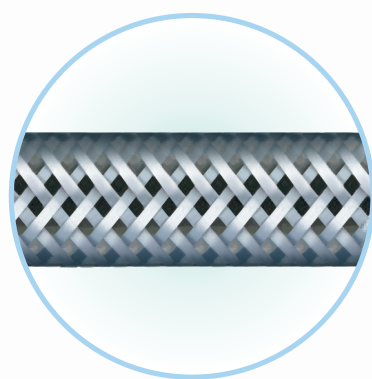

### Next Generation Hydrophilic Coating

Durable lubricity provides excellent crossability and trackability



### Stainless Steel Braided Mesh

Provides guidewire support for improved pushability

### Floppy Distal Segment

Distal 15 cm flexible length improves trackability through tortuous anatomy

### PTFE Inner Layer

Reduces friction and improves guidewire handling

### Gold Radiopaque Marker

For enhanced visibility that facilitates accurate placement

### Soft, Atraumatic Distal Tip

1.7 Fr entry profile, tapered distal tip for optimal guidewire handling and improved crossability

## PROVEN TECHNIQUES TO OPTIMIZE YOUR APPROACH

Access and cross simple to complex lesions using FINECROSS® M3 Coronary Micro-Guide Catheter with RUNTHROUGH® NS Coronary Guidewires


### SUPPORT for Side-branch or Acute Bifurcation Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Extra Floppy**

Coronary Guidewires

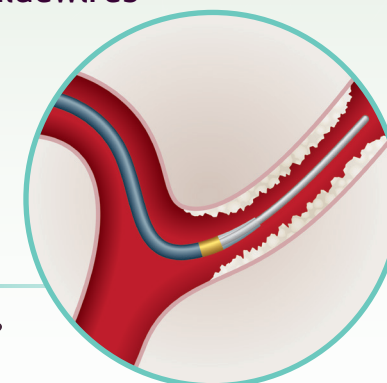


### TRACKABILITY for Distal or Torturous Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Izanai™**

Coronary Guidewires

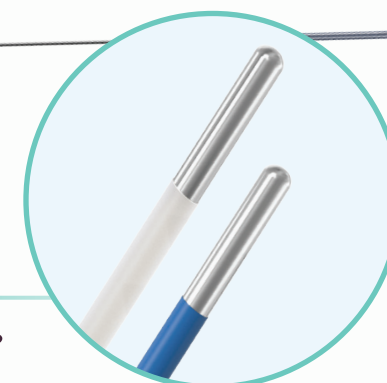

RUNTHROUGH NS Izanai is available in blue and white for simple identification of secondary wires


### CROSSABILITY for Tight Stenotic or Calcified Lesions

Use **FineCross® M3**

Coronary Micro-Guide Catheter

with **Runthrough® NS Hypercoat™**

Coronary Guidewires

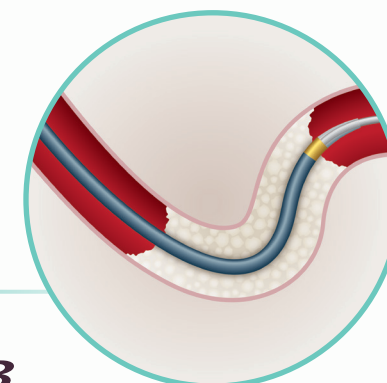

**TERUMO**  
INTERVENTIONAL  
SYSTEMS

# FineCross® M3

Coronary Micro-Guide Catheter

# Runthrough® NS

Coronary Guidewires

## ORDERING INFORMATION

### FINECROSS® M3 Coronary Micro-Guide Catheter

PRODUCT CODE	EFFECTIVE CATHETER LENGTH (cm)	DISTAL ENTRY TIP PROFILE (FR/MM)	DISTAL OUTER DIAMETER (Fr/mm)	PROXIMAL OUTER DIAMETER (Fr/mm)	DISTAL HYDROPHILIC COATING (CM)	DISTAL INNER DIAMETER (MM)	PROXIMAL INNER DIAMETER (IN/MM)	MAXIMUM GUIDEWIRE SIZE (IN/MM)
M3-1430	130	1.7/0.55	1.8/0.61	2.5/0.84	75	0.45	0.021/0.55	0.014/0.36
M3-1450	150	1.7/0.55	1.8/0.61	2.5/0.84	95	0.45	0.021/0.55	0.014/0.36

### RUNTHROUGH® NS Extra Floppy Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-1011	0.014	180	3	Straight/Shapeable
25-1013	0.014	300	3	Straight/Shapeable

### RUNTHROUGH® NS Izanai™ Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE / SHAFT COLOR
25-5011	0.014	180	3	Straight/Blue
25-5211	0.014	180	3	Straight/White
25-5213	0.014	300	3	Straight/White

### RUNTHROUGH® NS Hypercoat™ Coronary Guidewire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-3011	0.014	180	3	Straight/Shapeable
25-3013	0.014	300	3	Straight/Shapeable

Packed 5 per box. Shaping tool included.

### RUNTHROUGH® NS Extension Wire

PRODUCT CODE	DIAMETER (in)	LENGTH (cm)	TIP RADIOPACITY (cm)	TIP SHAPE
25-9092	0.014	150	—	—

Packed 1 per box. Does not include a shaping tool.

**FIND OUT MORE**  Phone: 800-888-3786  [terumo.com](https://www.terumo.com)

#### References:

1. Data on file.
2. Claims limited to tested microcatheters. (Volcano Valet & ASAHI Corsair)

**RX ONLY.** Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

©2023 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of their respective owners. PM-06949

