

Glidewire Advantage Track™

Peripheral Guidewire

LEVERAGE SUPERIOR PERFORMANCE

For your demanding CLI procedures

Available in 0.018" & 0.014"
diameter sizes



The addition of **GLIDEWIRE ADVANTAGE TRACK™** Guidewire to your PAD/CLI toolkit provides you with **superior performance^{1,a}** for your most challenging cases:

- **Superior Pushability:** High deformation resistance provides superior pushability to access and cross challenging above-the-knee (ATK) & below-the-knee (BTK) anatomy and lesions
- **Superior Trackability:** Provides superior ability to track and navigate anatomy of the vasculature
- **Best-in-Class Durability:** Provides superior shape retention for navigating small, tortuous vessels
- **Superior Lubricity:** TERUMO Glide Technology™ hydrophilic coating^b provides superior lubricity to support ease of progression through vessels
- **Enhanced Visibility & Selectivity:** Gold coil tips enhance visibility in conjunction with a 25 cm distal taper designed to minimize vessel trauma and improve selectivity of vessels

^a When compared to Abbott Hi-Torque Command 18 LT and Boston Scientific V-18

^b Referring to the distal portion of GLIDEWIRE ADVANTAGE TRACK™ wire which is GLIDEWIRE®

PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.


We are relentlessly seeking new ways to help you achieve **better outcomes for more patients.**



GLIDEWIRE ADVANTAGE TRACK™ Peripheral Guidewire

PRODUCT CODE	DIAMETER	TOTAL LENGTH	SPIRAL PTFE PROXIMAL COATING LENGTH	DISTAL GLIDEWIRE® LENGTH	FLEXIBLE TAPER LENGTH (TIP)	TIP SHAPE
GAT1418	0.014"	180 cm	155 cm	25 cm	1 cm	35° Angle
GAT1430	0.014"	300 cm	275 cm	25 cm	1 cm	35° Angle
GAT1818	0.018"	180 cm	155 cm	25 cm	1 cm	35° Angle
GAT1830	0.018"	300 cm	275 cm	25 cm	1 cm	35° Angle

Technical Specifications

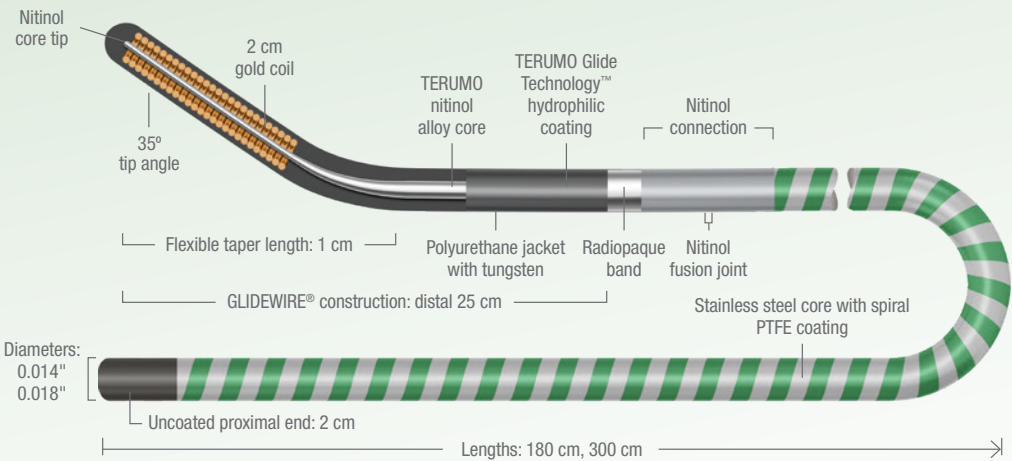

LENGTHS	DIAMETERS	FLEXIBLE TAPER LENGTH	TIP SHAPE	DISTAL MARKER	OTHER
180 cm 300 cm	0.014" 0.018"	1 cm	35° Angle	2 cm gold coil	25 cm GLIDEWIRE® distal end - Spiral PTFE-coated stainless steel core

FIND OUT MORE



US: 800.888.3786
CA: 833.883.7866



US: [terumo.com](https://www.terumo.com)
CA: [terumocanada.ca](https://www.terumocanada.ca)

Reference:
1. Data on file.

RX ONLY. The advertisement is directed to physicians only, and not to consumers. Refer to product labels and packaging insert for complete warnings, precautions, potential complications, and instructions for use.

©2024 Terumo Medical Corporation. All rights reserved. The third-party trademarks used herein are trademarks of their respective owners. All other brand names are trademarks or registered trademarks of Terumo. PM-08756

