

Glidewire®

Hydrophilic Coated Guidewire

MASTER YOUR CHOICE FOR PROCEDURE SUCCESS

The TERUMO GLIDEWIRE® product line that you know and trust offers an expanded range of tip shapes and wire configurations, making it possible for you to find the right tools for your challenging procedures.

1. Stiff GLIDEWIRE®

- Extra support and guidewire pushability designed to help straighten tortuous vasculature for better device delivery

2. 1 cm Taper GLIDEWIRE®


- Shorter taper designed to help cross stenotic lesions

3. Long Taper GLIDEWIRE®

- Long taper designed to both help maintain catheter tip shape and gain deeper vessel purchase to support guidewire advancement

4. J-Tip GLIDEWIRE®

- J-curve designed to provide atraumatic tip and avoid selection of vessel side branches

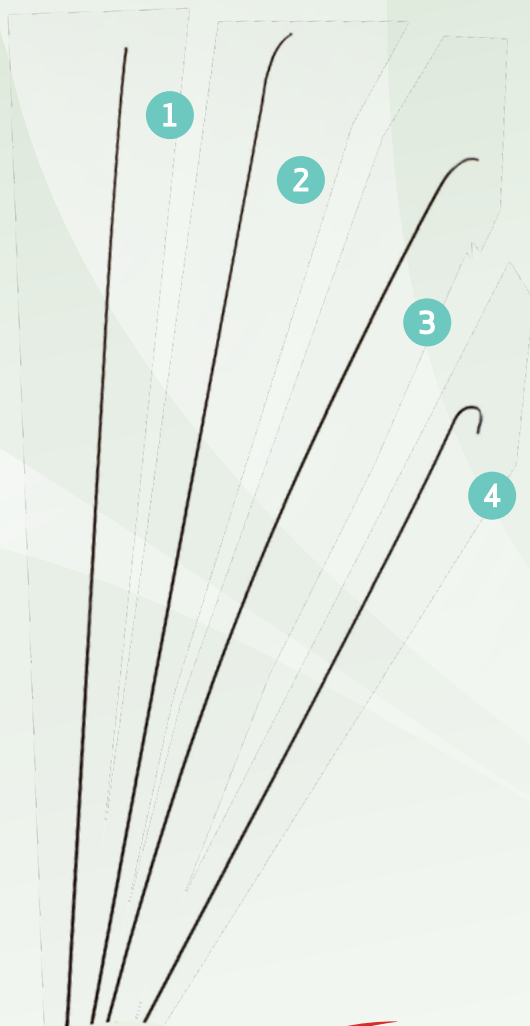

TERUMO
INTERVENTIONAL
SYSTEMS

PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you achieve **better outcomes for more patients**.

Consistent, reliable TERUMO GLIDEWIRE® performance designed to get you there

With TERUMO Glide Technology™ hydrophilic coating proven to have superior lubricity retention^{1-3*}, and a super-elastic nitinol alloy core for optimal performance and kink resistance⁴, the TERUMO GLIDEWIRE® offers reliable performance you know and trust.

	PRODUCT CODE	WIRE TYPE	DIAMETER	FLEXIBLE TIP LENGTH (TAPER)	TIP SHAPE	TOTAL LENGTH
Stiff Shaft GLIDEWIRE®	GS3505	Stiff	0.035"	3 cm	Angle	80 cm
	GS3506					150 cm
	GS3508					180 cm
	GS3509					260 cm
	GS3501				Straight	150 cm
	GS3503					180 cm
	GS3504		260 cm			
	GS3804		Angle		80 cm	
	GS3805				150 cm	
	GS3807				180 cm	
	GS3801	Straight			150 cm	
	GS3803		180 cm			
Long Taper GLIDEWIRE®	GS3507	Stiff	0.035"	5 cm	Angle	150 cm
	GS3510				Straight	150 cm
	GS3502		0.038"		Angle	
	GS3806				Straight	
	GR3507	Standard	0.035"	8 cm	Angle	
	GR3502				Straight	
1 cm Taper GLIDEWIRE®	GR3515	Standard	0.035"	1 cm	Angle	150 cm
	GR3516					180 cm
	GR3517					260 cm
	GR3518				Straight	150 cm
	GR3519					180 cm
	GR3520					260 cm
J-Tip GLIDEWIRE®	GR3510	Standard	0.035"	3 cm	3.0 mm J-tip	150 cm
	GR3810		0.038"			1.5 mm J-tip
	GR3525		0.035"		260 cm	
	GR3526					

FIND OUT MORE ☎ Phone: 800.888.3786 🌐 [terumo.com](https://www.terumo.com)

*When compared to ZIPwire™, HiWire®, AQUATRACK®, and Laureate®.

References:

- 1. Shah A, Lau C, Stavropoulos SW, et al. Comparison of physician-rated performance characteristics of hydrophilic-coated guidewires. *J Vasc Interv Radiol*. 2008 Mar;19(3):400-405.
- 2. Niazi K, Farooqui F, Devireddy C, et al. Comparison of hydrophilic guidewires used in endovascular procedures. *J Invasive Cardiol*. 2009 Aug;21(8):397-400.
- 3. Niazi K. Guidewire Satisfaction Study. Data on File 2012, Terumo Medical Corporation.
- 4. Data on file.

RX ONLY. Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

©2024 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of their respective owners. PM-07704

