



Pinnacle® Destination®

Guiding Sheath

GUIDE YOUR
PROCEDURES
THROUGH
**CONSISTENT
RELIABILITY**

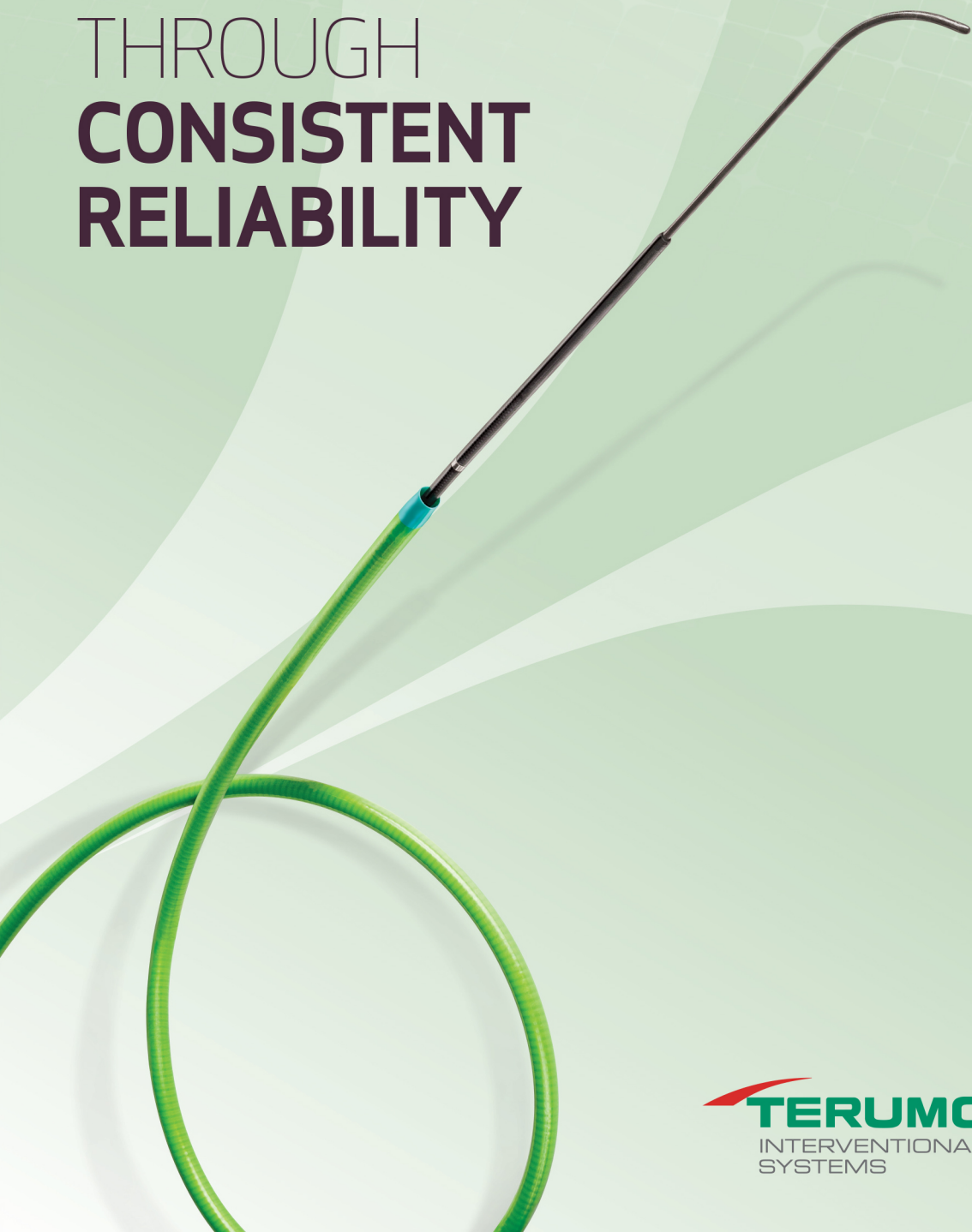

 **TERUMO**
INTERVENTIONAL
SYSTEMS

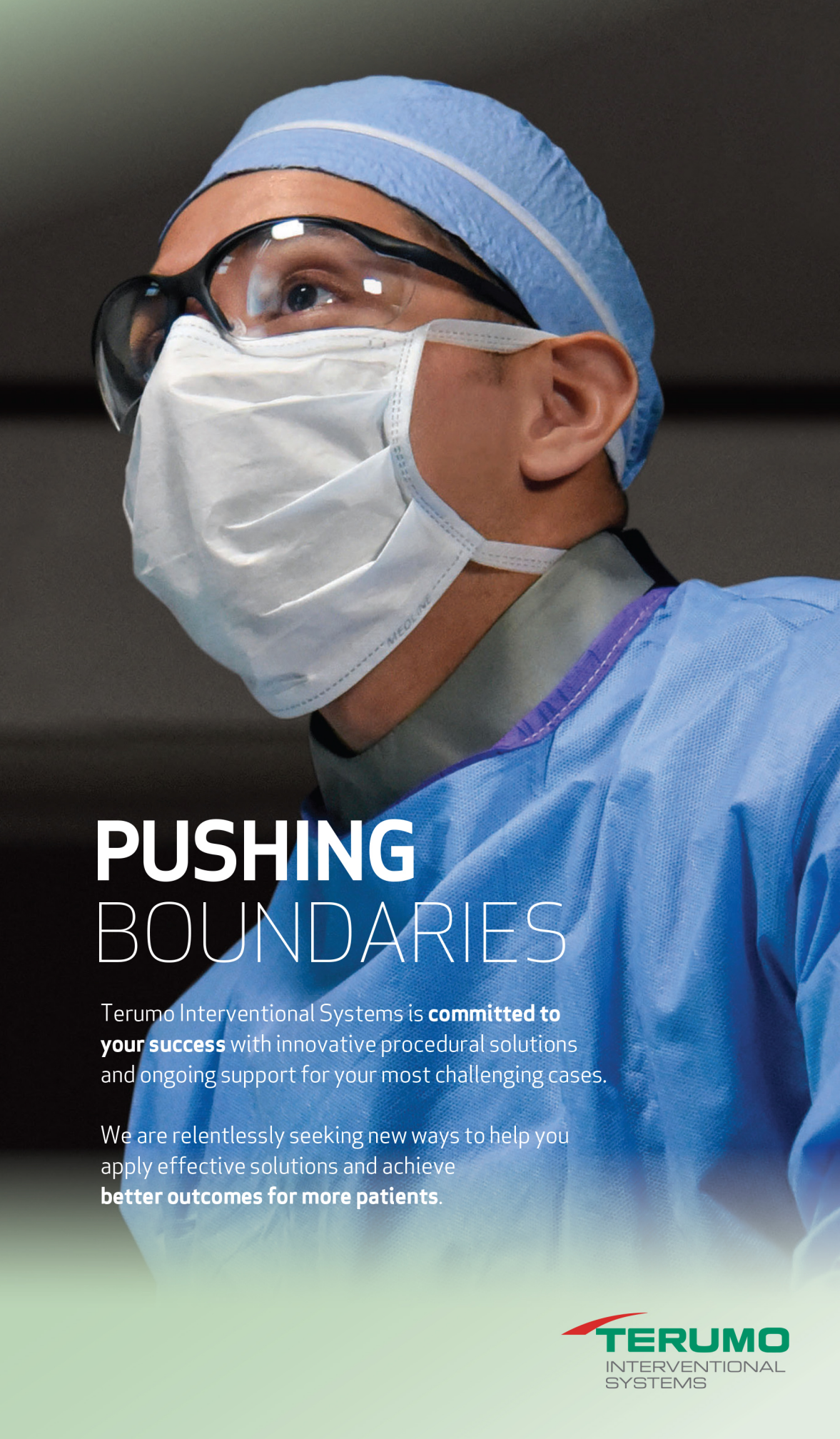

PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you apply effective solutions and achieve **better outcomes for more patients.**

Pinnacle® Destination®

Guiding Sheath

ALLOWS SMOOTH TRANSITION AND EASY PENETRATION

Atraumatic tip¹


- Minimizes potential for vessel damage

Smooth transitions¹

- Guidewire-to-dilator and dilator-to-sheath

Provides lesion access¹

- Exceptional trackability, even in the most challenging anatomy

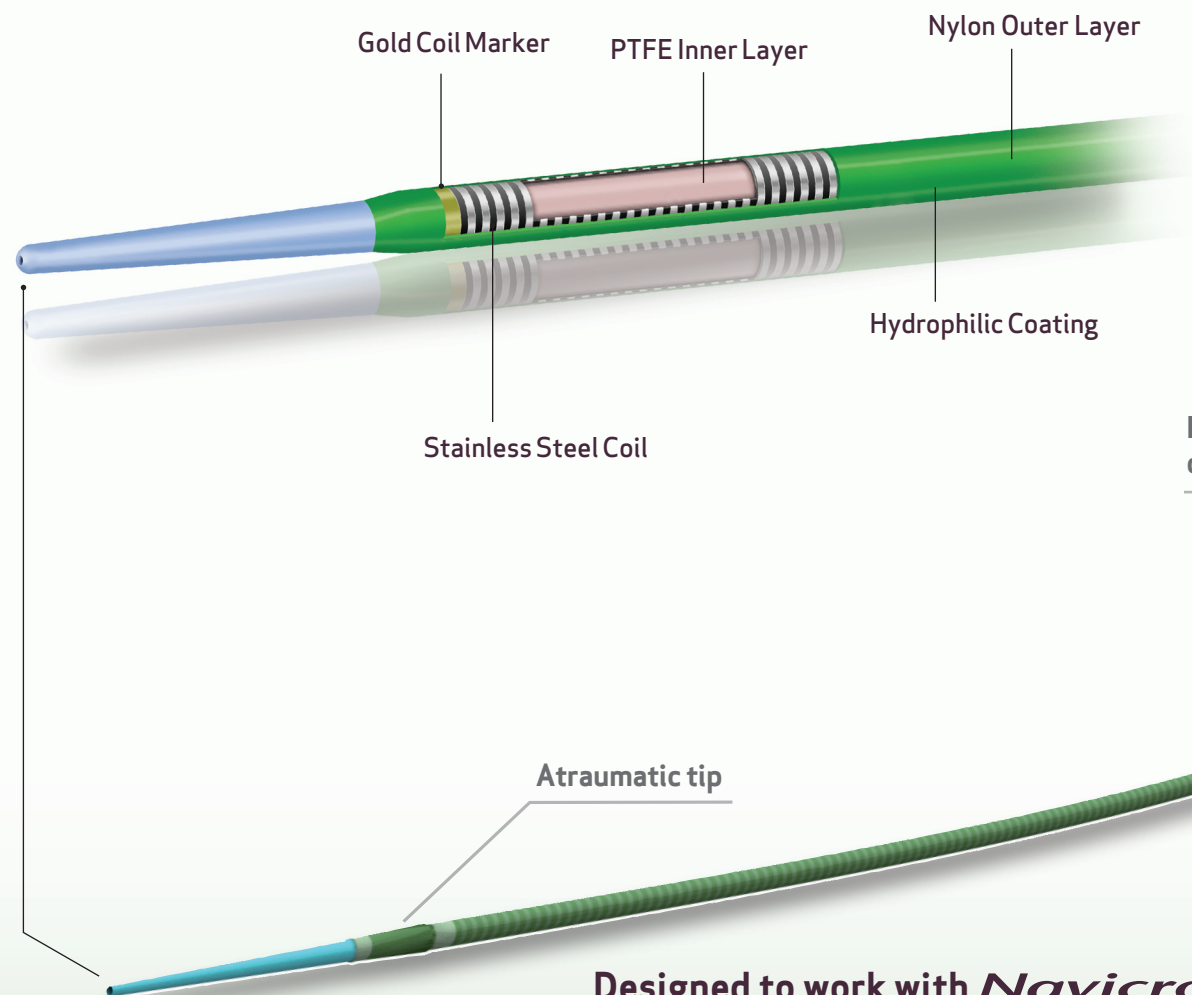

Designed to work with *Navicross*® and *Glidewire Advantage*®

Support Catheters

Peripheral Guidewires

SHAPES AND SIZES TO MEET YOUR DIAGNOSTIC CHALLENGES

Maintains shape and positioning¹

- Outstanding durability, even in the longest procedures

Enhances visualization for precise positioning¹

- Triple radiopacity (sheath, dilator, and gold coil marker)

Resists kinking¹

- Stainless steel coil for consistent reliability

Provides smooth movement¹

- Coil reinforced tubing and a PTFE inner layer that minimizes friction

Detachable valve options

Nylon outer layer reinforced with stainless steel coil

Hydrophilic coating on distal end

Cross Cut™ valve


Tuohy-Borst valve

Accommodates anatomical variations in vessel access

- Straight
- Hockey Stick
- Multi-Purpose
- RDC (Renal Double Curve)
- LIMA (Left Internal Mammary Artery)

Pinnacle® Destination®

Guiding Sheath

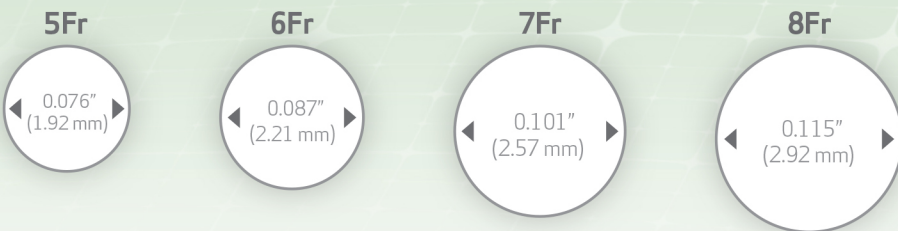

PINNACLE® DESTINATION® Guiding Sheaths

Packed 1 per box and includes dilator. All dilators are 0.038" wire compatible.

French Size	Product Code	Length (cm)	Hydrophilic Coating Length (cm)	Valve Type	Curve Style	Curve Shape
5Fr	54-54501	45	35	Cross-Cut	Straight	
	54-54502	45	5	Cross-Cut	Hockey Stick	
	54-54503	45	5	Cross-Cut	Multi-Purpose	
	54-54504	45	5	Cross-Cut	RDC	
6Fr	54-64501	45	35	Cross-Cut	Straight	
	54-64506	45	35	Tuohy-Borst	Straight	
	RSR01	45	5	Cross-Cut	Straight	
	RSR02	45	5	Cross-Cut	Hockey Stick	
	RSR03	45	5	Cross-Cut	Multi-Purpose	
	RSR07	45	5	Tuohy-Borst	Straight	
	RSR08	45	5	Tuohy-Borst	Hockey Stick	
	RSR09	45	5	Tuohy-Borst	Multi-Purpose	
	RSR13	45	5	Cross-Cut	RDC	
	RSR14	45	5	Cross-Cut	LIMA	
	RSR17	45	5	Tuohy-Borst	RDC	
	RSR18	45	5	Tuohy-Borst	LIMA	
	RSP01	65	35	Cross-Cut	Straight	
	RSP03	65	35	Tuohy-Borst	Straight	
	RSC01	90	15	Tuohy-Borst	Straight	
	RSC03	90	15	Tuohy-Borst	Multi-Purpose	
7Fr	54-74501	45	35	Cross-Cut	Straight	
	54-74506	45	35	Tuohy-Borst	Straight	
	RSR04	45	5	Cross-Cut	Straight	
	RSR05	45	5	Cross-Cut	Hockey Stick	
	RSR06	45	5	Cross-Cut	Multi-Purpose	
	RSR10	45	5	Tuohy-Borst	Straight	
	RSR11	45	5	Tuohy-Borst	Hockey Stick	
	RSR12	45	5	Tuohy-Borst	Multi-Purpose	
	RSR15	45	5	Cross-Cut	RDC	
	RSR16	45	5	Cross-Cut	LIMA	
	RSR19	45	5	Tuohy-Borst	RDC	
	RSR20	45	5	Tuohy-Borst	LIMA	
	RSP02	65	35	Cross-Cut	Straight	
	RSP04	65	35	Tuohy-Borst	Straight	
	RSC02	90	15	Tuohy-Borst	Straight	
	RSC04	90	15	Tuohy-Borst	Multi-Purpose	
8Fr	54-84501	45	35	Cross-Cut	Straight	
	54-84506	45	35	Tuohy-Borst	Straight	
	54-86501	65	35	Cross-Cut	Straight	
	54-86506	65	35	Tuohy-Borst	Straight	
	54-89001	90	60	Cross-Cut	Straight	
	54-89006	90	60	Tuohy-Borst	Straight	

FIND OUT MORE | Phone: 800.888.3786 [terumo.com](https://www.terumo.com)

RX ONLY. Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

Reference:
1. Data on File.

©2023 Terumo Medical Corporation. All rights reserved.
All brand names are trademarks or registered trademarks of Terumo. PM-06999

